

Reorder Nos. 10750-10752

 **Manufactured for:**
Dynarex Corporation
10 Glenshaw Street
Orangeburg, NY 10962
USA • www.dynarex.com

Made in China



Standard Sling




SYMBOL GLOSSARY
For an explanation of symbols
used in Dynarex packaging, visit
dynarex.com/symbols.php

WARNINGS

- Read the instructions before using the sling.
- Before use consult with a trained professional to ensure the proper sling is used. Improper use of the sling may result in personal injury.
- Carefully inspect the sling for wear before each use. If the sling shows signs of tearing, fraying or other wear, discard the sling immediately. Damaged slings are not safe for use and may result in serious personal injury to the user.
- Do not exceed maximum weight capacity. Never exceed the lowest weight capacity rating of either the patient lift or patient sling.
- This sling should be used with Dynarex brand lift, and chains or straps only.
- Use only with chains or straps with “S” shaped hooks.

ASSEMBLY INSTRUCTIONS

1. Locate the “O” shaped ring on the chain.
2. Connect the “S” shaped hook on the shorter chain to the back panel of the standard sling.
3. Connect the “S” shaped hook on the longer chain to the front panel of the standard sling.
4. Connect the “O” shaped ring to the Lift.

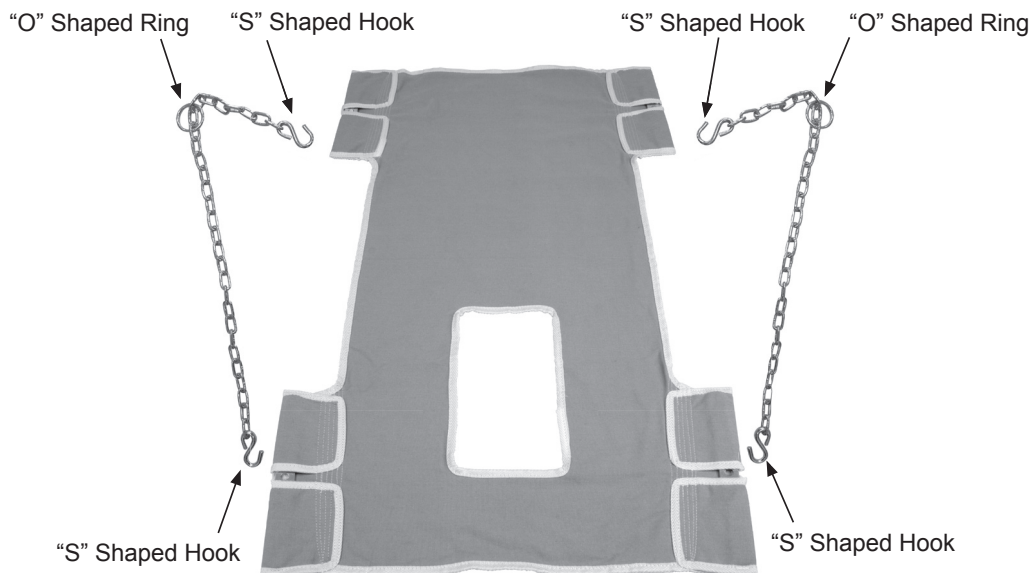

SIX MONTH LIMITED WARRANTY

Your Dynarex Product is warranted to be free of defects in materials and workmanship for six (6) months from original date of purchase, for the original purchaser. This item was built to exacting standards and carefully inspected prior to shipment. This Warranty is an expression of our confidence in the materials and workmanship of our products and our assurance to the consumer of years of dependable service.

The Warranty shall not apply under the following conditions:

- Problems arising from normal wear
- Problems arising from failure to adhere to the product instructions
- Problems arising from misuse, negligence, accident or improper operation, maintenance or storage
- Problems arising from modifications or unauthorized repairs, parts or attachments

- Products where the serial number has been removed or defaced
- Problems with non-durable components, such as rubber accessories, casters, and grips, which are subject to normal wear and need periodic replacement

Dynarex shall not be liable for any consequential or incidental damages whatsoever. Dynarex shall repair or replace defective products at its option. The foregoing warranty is exclusive and in lieu of other express warranties, if any, including the implied warranties of merchantability and fitness of a particular purpose. The remedy for any violation of the implied warranty shall be limited to repair or replacement of the defective product pursuant to the terms contained herein.

If you have a question about your Dynarex device or this warranty, please contact an authorized Dynarex dealer.