




Clearview® RSV

Rapid technology provides direct detection of the Respiratory Syncytial Virus (RSV)

The Clearview® RSV test is patented, innovative technology that provides highly sensitive and specific tests for Respiratory Syncytial Virus.

- **Rapid Response:** 4 easy steps, result in 15 minutes.
- **Simplicity:** CLIA waived, easy-to-use test requires no training.
- **Reliability:** Built-in control.
- **Flexibility:** Multiple sample types — nasopharyngeal swabs or aspirates.

The Clearview® RSV test is designed to aid in the rapid diagnosis in children less than 6 years and adults over age 60.

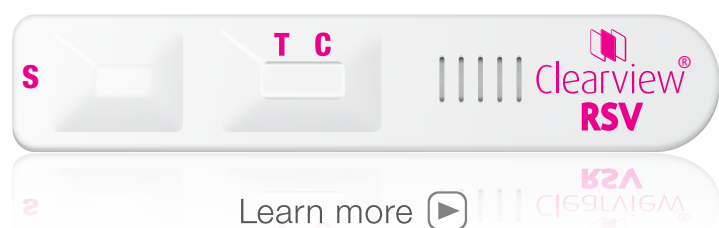


Clearview® RSV Comparison Grid

Product Name	Clearview® RSV	Quidel QuickVue® RSV	Alere BinaxNOW® RSV
Format	Cassette	Dipstick	Patented card format
CLIA Waived	YES	YES	YES
Sample Types	NP swab NP aspirate	NP swab NP / Nasal wash NP aspirate	NP swab Nasal wash
Reagents	1	1	0
Timed Steps	1	2	1
Tests per Kit	20	20	22
Shelf Life ¹	12 Months	24 Months	24 Months
Storage	Room Temperature	Room Temperature	Room Temperature
QC Requirements ²	Each new lot or shipment	Every new operator, each new lot or every shipment (provided that each different lot in the shipment is tested)	Once per shipment

Procedure

1

Gently agitate the nasopharyngeal swab or aspirate sample. For swab sample, roll the swab against the inside of the transport tube or specimen container as you remove it and discard the used swab.

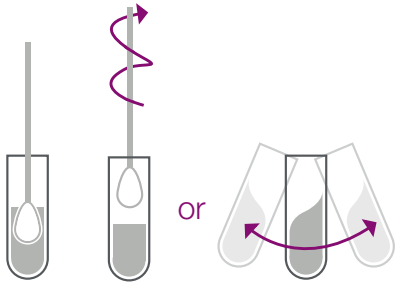


2

Fill the disposable pipette to the second notch (160ul) with the nasopharyngeal swab solution or nasal aspirate and add entire contents of the Extraction Tube.




3

Add four (4) free falling drops of R1 to the Extraction Tube while holding the bottle vertically over the Tube and gently swirl the Extraction Tube to mix the sample and R1.

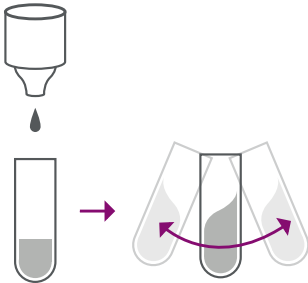

4

Draw entire content of the Extraction Tube up into a clean Disposable Pipette and squeeze the entire contents onto the sample well of the Test Device.

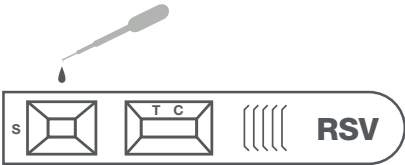


5


The test results should be read between fifteen (15) and twenty-five (25) minutes after applying the extracted sample to the Test Device.





6

Results:

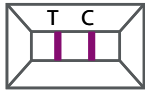
Positive

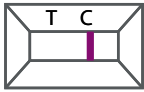
Negative

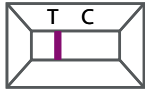
Invalid

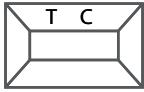
Invalid



Call **1.877.441.7440** or visit **alere.com** today

Clearview® RSV
20 Test Kit
Order No. **135060**
CIT: 87807QW

CLIA Waived

1. From date of manufacture 2. See package inserts for details
© 2013 Alere. All rights reserved. The Alere Logo, Alere, BinaxNOW and Clearview are trademarks of the Alere group of companies. All other trademarks referenced are trademarks of their respective owners. 12000289E-02 12/13