

TECHNICAL DATA SHEET



FEATURES & BENEFITS

- Thoughtful designs offer ease of movement, as well as secure closure and easy removal.
- Special features such as full back closure, string belt, knitted cuffs, and tie-on neck closure.
- Made from Spunbond Polypropylene, Polyethylene coated.
- Premium materials provide superior protection, reliability and durability.
- Lightweight, soft, supple, fabrics keep wearers cool and comfortable.

RECOMMENDED USE

- Procedure Isolation Gown designed to protect health professionals and patients from cross-contamination. Single use and disposable gowns. Latex-free.

PRODUCT DETAILS

Product Name:	Medicom SafeBasics™ Fluid Impervious Gown
Product Type:	Single-use, non-sterile
Color:	White
Sizes:	See Appendix 1
Quantity per Bag:	10
Country of Origin:	China

DONNING INSTRUCTIONS

- Perform hand hygiene before removing item from package.
- Select appropriate garment for intended use (length, coverage, fluid resistance, etc.)
- Fasten correctly: tied at neck and waist to ensure comfortable fit, adequate coverage and freedom of movement.
- Protective apparel should not be worn in restrooms, personal offices or break rooms.

DOFFING INSTRUCTIONS

- Remove apparel immediately after use. Do not re-use.
- Perform proper hand hygiene after removing apparel to prevent cross-contamination.
- Turn garment inside out when removing to keep infectious material inside.
- Designate a safe receptacle in which to dispose of protective apparel after use.

CERTIFICATION & STANDARDS


ISO 13485 certified manufacturer.
Class 1 Medical Device.
Health Canada registered medical device manufacturer.

STORAGE CONDITIONS

These gowns should be disposed of after each use.
Store under ambient temperature and dry condition.



GOWN SPECIFICATIONS			
CHARACTERISTIC		SPECIFICATION	
Material Composition		Spunbond Polypropylene, Polyethylene Coated	
Description		Gown made of non-woven fabric with full back closure. Features: String belt, knitted cuffs, and tie-on neck closure.	
Part of Gown		Dimension (cm)	
		Large	X-Large
A	Gown Length	110	115
B	Underarm Width	71	77,5
	Open Gown Width	142	155
C	Sleeve Length	54	58
D	Shoulder Height*	28	32
E	Shoulder Width	23	26
F	Knitted Cuff Width	6	6
G	Belt Length	180	180
	Belt Position	24 cm below the arm, sewn on a distance of 4 - 6 cm	24 cm below the arm, sewn on a distance of 4 - 6 cm
I	Collar Depth	9	9
J	Collar Width	23	23
L	Sleeve Cuff Width (Stretched)	15	15
	Knitted Cuff Length	6	6
	Neck-tie Length	35	35

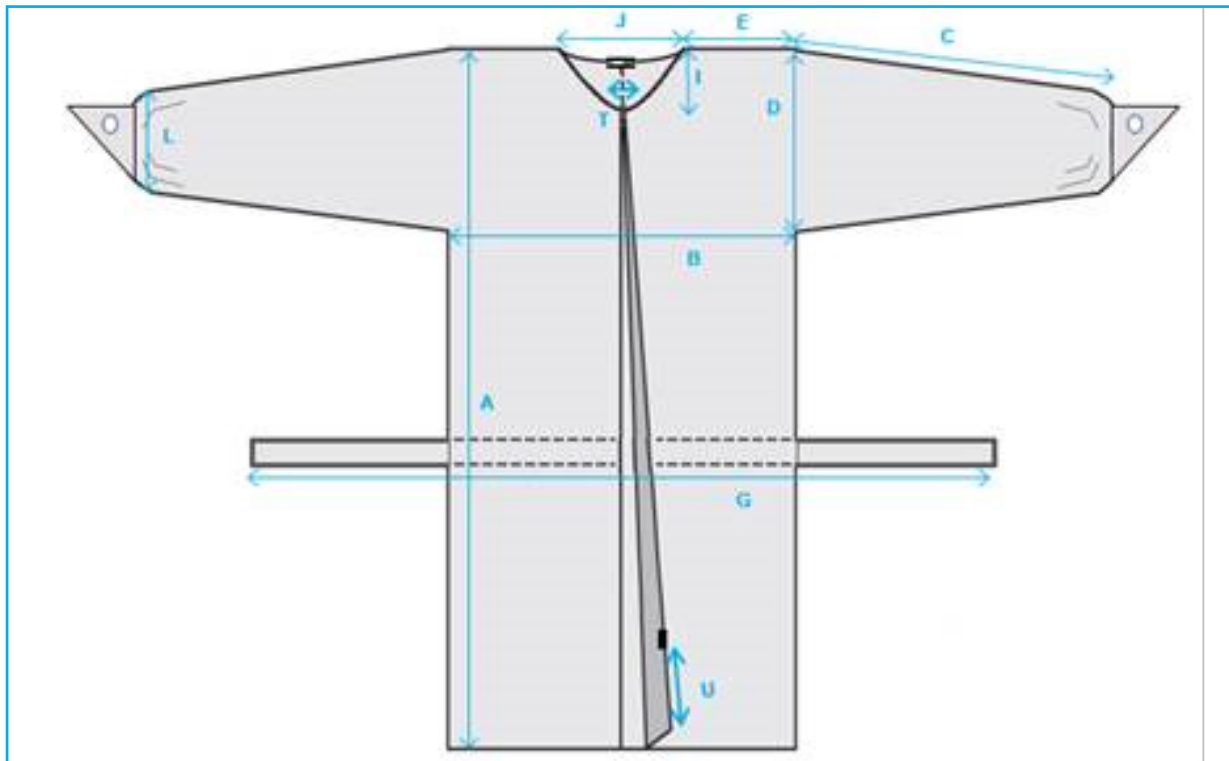

Appendix 1 (List of items)

Item Code	Description	Packaging Information	Size	Color	Specific Attribute	Treatment	Marketed in USA as	Marketed in Canada as
8022	Medicom SafeBasics Fluid Impervious Gown	10 Gowns/Bag, 5 Bags/Case, 50 Gowns/Case	Regular	White	Belt on the right side	Non-Sterile	Medical Device Class 1	Medical Device Class 1
8024	Medicom SafeBasics Fluid Impervious Gown	10 Gowns/Bag, 5 Bags/Case, 50 Gowns/Case	Extra Large	White	Belt on the right side	Non-Sterile	Medical Device Class 1	Medical Device Class 1

AMD Medicom Inc. | 2555 Chemin de l'Aviation, Montreal, Quebec, Canada H9P 2Z2