



SYMBOL GLOSSARY
For an explanation of symbols
used in Dynarex packaging, visit
dynarex.com/symbols.php

Overbed Table

White U Base



Reorder No. 10453
Reorder No. 10454
Reorder No. 10455

Important Warnings

Please follow all instructions and warnings to avoid hazardous situations that can result in property damage and/or patient/resident/staff/caregiver injury and/or death. Inspect the product for damage and loose hardware before each use.

DO NOT use the product if it appears to be damaged or unstable.

DO NOT modify the product.

Intended Use

This product is designed for use as a movable patient/resident table top surface for reading, eating, or keeping small personal items such as glasses, telephone, or books within easy reach. This product is designed for indoor use **ONLY** within close proximity of skilled caregivers. It is not designed to hold electrical items (such as televisions, computer monitors, coffee pots, etc.) and should not be used near open flames and excessive heat.

Maximum weight limitation for this product is 50 lb. (22.7 kg.); weight should be evenly distributed over the table top surface while at the lowest usable height to help maintain a stable condition. **DO NOT** overload the product as that could result in breakage, as well as damage to property and/or serious injury. Use **ONLY** as directed for its intended use.

This product should **NOT** be used: as a sleep surface; as a walker or mobility aid to assist a person in walking, sitting down or standing up; to prevent or brace a person's fall; or to help transport persons or items.

Limited Lifetime Warranty


Your Dynarex product is warranted to be free of defects in materials and workmanship for the lifetime of the product for the original consumer purchaser. This item was built to exacting standards and carefully inspected prior to shipment. This Limited Lifetime Warranty is an expression of our confidence in the materials and workmanship of our products and our assurance to the consumer of years of dependable service.

The Warranty shall not apply under the following conditions: problems arising from normal wear; problems arising from failure to adhere to the product instructions; problems arising from misuse, negligence, accident or improper operation, maintenance or storage; problems arising from modifications or unauthorized repairs, parts or attachments; products where the serial number has been removed or defaced; problems with non-durable components, such as rubber accessories, casters, and grips, which are subject to normal wear and need periodic replacement. Dynarex shall not be liable for any consequential or incidental damages whatsoever. Dynarex shall repair or replace defective products at its option. The foregoing warranty is exclusive and in lieu of other express warranties, if any, including the implied warranties of merchantability and fitness of a particular purpose. The remedy for any violation of the implied warranty shall be limited to repair or replacement of the defective product pursuant to the terms contained herein.

Assembly Instructions


- 1) Remove contents from packaging and lay Table Top on a flat surface with mounting holes facing upward. (see Figure 1)

Figure 1

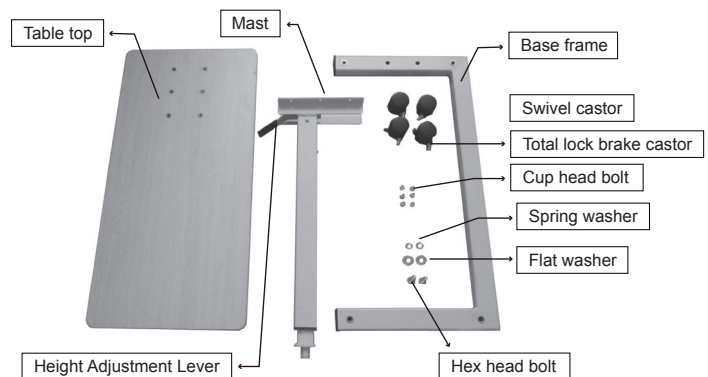

- 2) Secure the Mast to the Table Top with 6 Cup head bolts. Tighten securely.
- 3) Install the total lock brake castors diagonally. Insert the caster stem into the bolt hole in the base frame. Fasten the caster stem by tightening the bolt at the base of the caster. (see Figure 2)

Figure 2

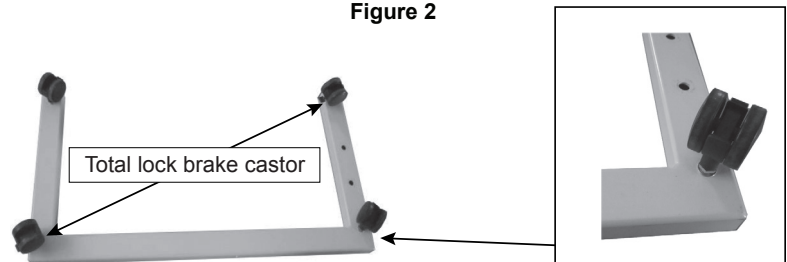

- 4) Secure the Base to the Mast by inserting a spring washer, flat washer, and Hex head bolt through the base into each bolt hole at the bottom of the mast and securing tightly. (see Figure 1)
- 5) Turn the table upright and stand on casters.
- 6) Remove black transport screw located on the Mast just below table top.

Operating Instructions

Adjusting Table Height:

- 1) Gently squeeze the Height Adjustment Lever up towards the table top. Adjust table top to desired height.
- 2) Check that Height Adjustment Lever located underneath the table top is properly engaged prior to use.