500HB Installation Instructions

PELSTAR, LLC 9500 West 55th St. McCook, IL 60525-7110 USA
© Pelstar, LLC 2020

www.homscales.com

weigh**easier**®

Thank you for purchasing this Health o meter[®] Professional product. To prevent patient/caregiver injury or damage to your scale, please follow the instructions in this user manual very carefully.


- The handlebars are intended for patient use only. Do not use the handlebars to transport the scale or use the scale as a handcart to transport items.
- To prevent patient injury, the patient must be attended throughout the entire weighing event. Caregivers should ensure that the patient is stable and provide support as needed when getting onto and off of the scale.
- Do not allow the patient to use the handlebars to pull themselves onto the scale.
- Excessive movement of the handlebars may affect the weight lock function.
- Do not allow patients to use excessive force to push or pull on the handlebars during the weighing process.

In no event whatsoever shall Pelstar[®], LLC be liable for damages or injuries whatsoever arising out of or connected with the assembly, use, or misuse of its products.

Tools Included: ¼" Drill Bit, Allen Wrench

Tools Required: Power Drill, Utility Knife

Parts List:

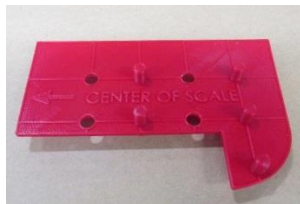Drill Guide

Alcohol Wipe

¼" Drill Bit

Allen Wrench

1/4-20 X 2 1/4
Button Head
Screws, 8 Each

1/4-20 X 1/2
Button Head
Screws, 4 Each

500HB Mount
Bracket, 2 Each

500HB SPACER

500HB Base

500HB Handlebar

Mounting Handlebars on New, Unassembled Scales

Refer to the scale's user manual for assembly instructions. In step 8 of the scale assembly, do not completely insert all 4 screws. Partially insert only the two screws closest to the scale column, leaving $\frac{3}{4}$ " of the screws exposed. Then skip to step 5 on the following page.

Mounting Handlebars on Assembled Scales


1. Gently place the scale on a flat, level surface on its side.


2. Remove the 2 Rear Screws (closest to the wheels) that secure the Scale Column to the Base.


3. Loosen but **DO NOT REMOVE** the 2 Forward Screws (closest to the scale base). Loosen Screws approximately $\frac{3}{4}$ ". *Note, the Load Cell Cable does not need to be removed.*


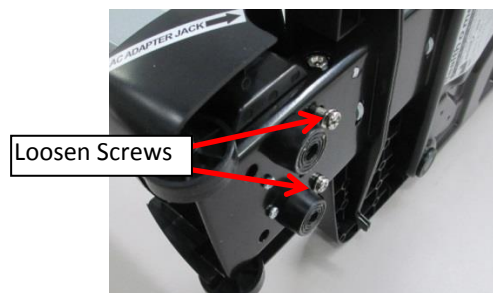
4. Loosen Screws approximately $\frac{3}{4}$ ".

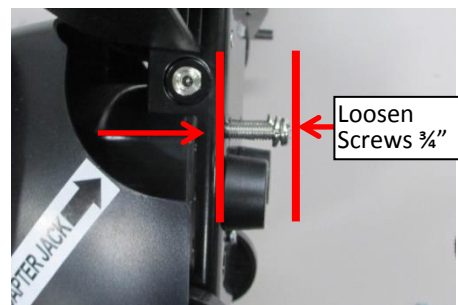
5. Tilt the front edge of the Base toward the top of the Scale opening a gap at the rear of the Scale. **Caution: Ensure the load cell cable does not disconnect from the load cell port.**

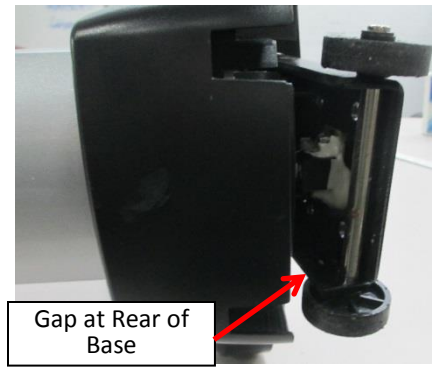6. Through the opening in the gap, slide the Spacer into the Gap. Position the Spacer so that the open end of the Spacer slides around the Load Cell Cable.

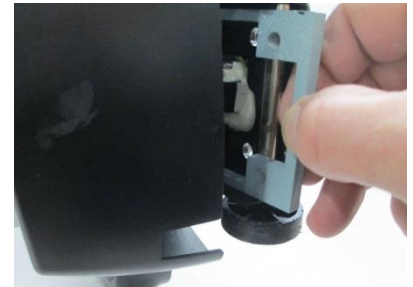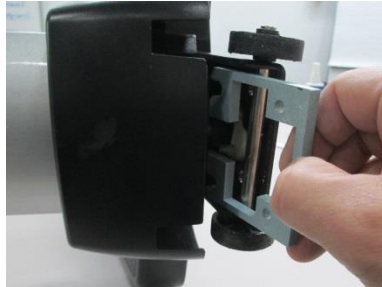
7. Slide Spacer into Base and position Spacer so that Mount Holes in Spacer align with Holes in Base.

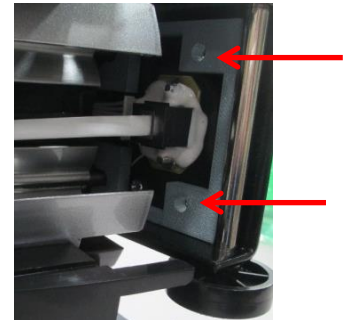
8. Insert the two Screws removed in step 3 into the Mount Holes on bottom of Base. Then tighten all 4 Screws securely.

Note: If assembling a new scale, refer back to the scale's user manual to complete the scale assembly before proceeding with mounting the handlebars.

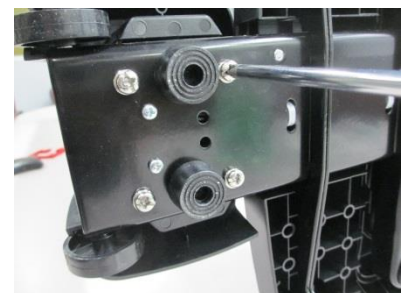
9. Place the drill guide onto the bottom of the scale platform as shown. Position drill guide so that the arrow on the guide points to center of the scale. Align edge of drill guide with edge of platform.

Align edge of drill guide with edge of scale platform

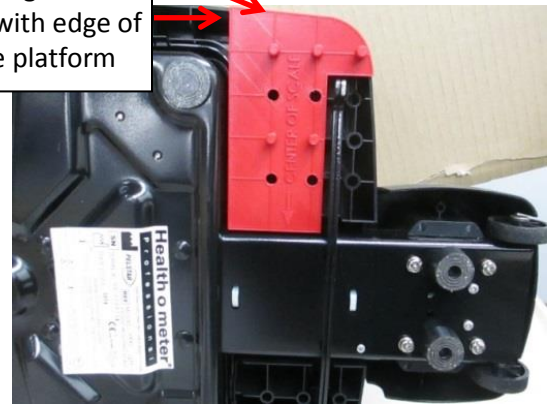
10. There are posts on both sides of the drill guide. These posts must slide into holes in platform.


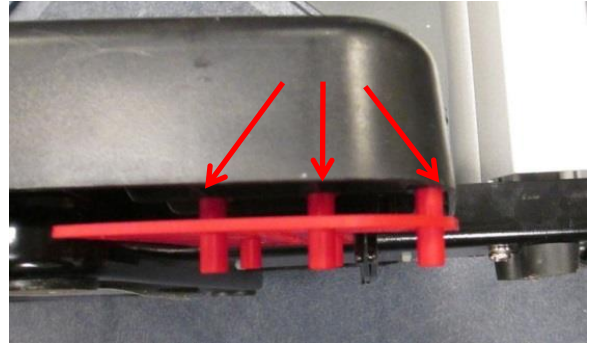
11. Drill guide must be positioned so that posts slide **fully** into holes on bottom of scale platform.

12. With drill guide properly positioned on bottom of scale platform, use a power drill and the provided $\frac{1}{4}$ " drill bit to drill 4 holes through the bottom of the scale platform.

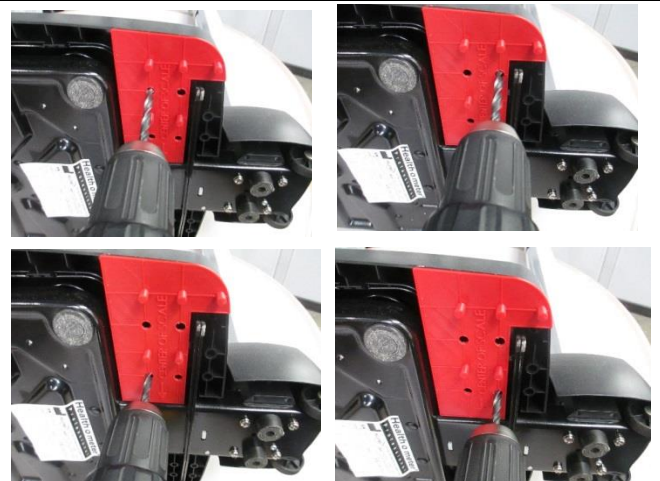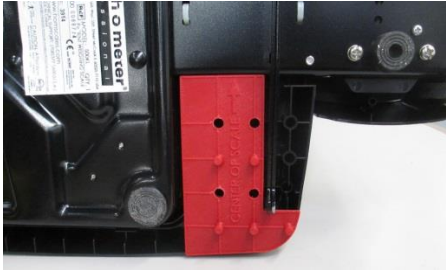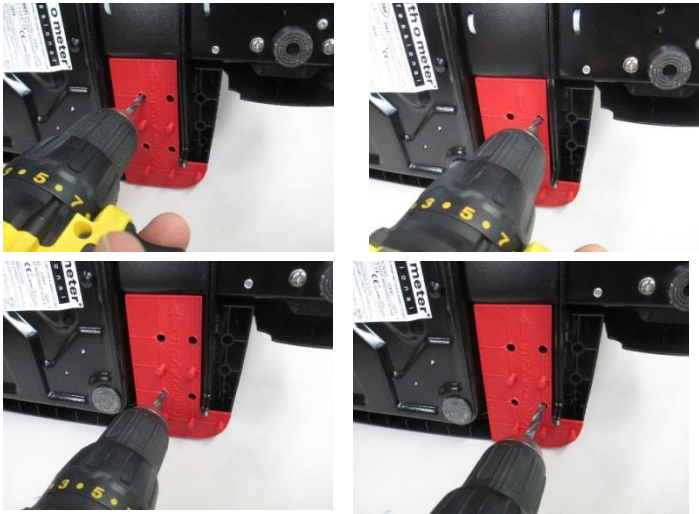13. Flip the drill guide (drill guide is used on both sides of the scale platform) and position it on the other side of the bottom of the platform.	
14. With drill guide properly positioned on bottom of scale platform, use a power drill and the provided ¼" drill bit to drill 4 holes through the other of the side bottom of the scale platform.	
15. Remove and discard drill guide. Carefully place the scale upright on the floor. Platform will now have 8 holes for mounting the handlebars.	
16. Use a utility knife or similar tool to remove plastic burrs from the top of the scale platform.	
17. Using the provided alcohol wipe, clean the entire top surface of the scale platform.	

18. Position the base onto the scale platform by aligning the holes in the plate with the holes in the drilled in the platform. Note: There is a strip of double-sided tape on the bottom of the handlebar base. This will be removed in step 22.

19. Insert two of the provided 1/4-20 X 2 1/4" screws into the plate on the handlebars and into the scale platform, one screw on the outside of both sides of the handlebar channel.

20. Locate the strip of double-sided tape on the bottom of the plate.

21. Slightly tilt the base up to allow for the removal of the tape backing strip. Peel off the backing strip under the plate. Lower the plate/handlebars onto the scale platform and apply light pressure for 30 seconds.

22. Insert the remaining six button head screws through the base and into the scale platform. Do not completely tighten.

23. Access the bottom of the platform by placing the scale carefully on its side or up on a table. Obtain one of the silver mounting brackets and position it so that the flat side of the bracket faces the bottom of the platform.

24. Position the bracket as shown and thread the screw into the threaded nut on the bracket.

25. Repeat the process with the second mounting bracket. Then tighten all 8 screws completely.

26. **Note:** When mounting the brackets underneath the platform, ensure that the silver brackets are not touching the black stabilizer bar. If necessary, loosen the screws to adjust, then retighten.

27. Obtain the handlebars and place around the scale pillar as shown.

28. Align the bottom of the handlebar with the channels on the base. Slide the handlebar onto the channel until the bottom of the handlebar rests on the base. Ensure the holes in the handlebar align with the holes in the channel.

29. Insert two of the ¼-20 x ½" Button Head Screws into the back of the handlebar and into the channel. Use the included Allen wrench to fully tighten and secure the screws.

30. Insert the remaining two ¼-20 x ½" Button Head Screws into the back of the other side of the handlebar and into the channel. Fully tighten to secure.

Assembly is complete.

WARRANTY

Limited Warranty

What does the Warranty Cover?

This Health o meter® Professional product is warranted from date of purchase against defects of materials or in workmanship for a period of two (2) years. If product fails to function properly, return the product, freight prepaid and properly packed to Pelstar, LLC (see "To Get Warranty Service", below, for instructions). If the manufacturer determines that a defect of material or in workmanship exists, the customer's sole remedy will be replacement of the product at no charge. Replacement will be made with a new or remanufactured product or component. If the product is no longer available, replacement may be made with a similar product of equal or greater value. All replaced parts are covered only for the original warranty period.

Who is Covered?

The original purchaser of the product must have proof of purchase to receive warranty service. Please save your invoice or receipt. Pelstar dealers or retail stores selling Pelstar products do not have the right to alter, or modify or any way change the terms and conditions of this warranty.

What is Excluded?

Your warranty does not cover normal wear of parts or damage resulting from any of the following: negligent use or misuse of the product, use on improper voltage or current, use contrary to the operating instructions, abuse including tampering, damage in transit, or unauthorized repair or alternations. Further, the warranty does not cover natural disasters, such as fire, flood, hurricanes and tornadoes. This warranty gives you specific legal rights, and you may also have other rights that vary from country to country, state to state, province to province or jurisdiction to jurisdiction.

To get Warranty Service make sure you keep your sales receipt or document showing proof of purchase. Call (+1) 800-638-3722 or (+1) 708-377-0600 to receive a return authorization (RA) number, which must be included on the return label. Attach your proof of purchase to your defective product along with your name, address, daytime telephone number and description of the problem. Carefully package the product and send with shipping and insurance prepaid to:

Pelstar, LLC

Attention R/A# _____

Return Department

9500 West 55th Street

McCook, IL 60525

PELSTAR, LLC

9500 West 55th St. McCook, IL 60525-7110 USA

1-800-638-3722 or 1-708-377-0600

PLEASE REGISTER YOUR PRODUCT FOR WARRANTY COVERAGE AT:

www.homscales.com

Health o meter® is a registered trademark of Sunbeam Products, Inc. used under license. Health o meter® Professional products are manufactured, designed, and owned by Pelstar, LLC. We reserve the right to improve, enhance, or modify Health o meter® Professional product features or specifications without notice.

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