

GRIPPER PLUS® Safety Needle

Needle Selection

We offer the GRIPPER PLUS® safety needle, in addition to a full line of GRIPPER® straight and 90° bent PORT-A-CATH® needles. PORT-A-CATH® needles are used to access PORT-A-CATH®, PORT-A-CATH® II, ProPort®, and P.A.S. PORT® implantable access systems.

LENGTH SELECTION

Needle length selection should only be made after an assessment of the patient's portal. The GRIPPER PLUS® safety needle length is measured from the plastic base to the tip of the needle. The foam is included in the measurement length.


NEEDLE LENGTH IS DETERMINED BY:

- The needle should be long enough to reach the reservoir floor of the portal when inserted
- The foam cushion of the GRIPPER PLUS® safety needle should lie flush against the skin surface; neither the skin nor the cushion should be compressed under the base of the needle. The safety arm should be resting flat against the base of the GRIPPER PLUS® safety needle.

GAUGE SELECTION

- GRIPPER PLUS® safety needles are available in 19, 20, and 22 gauge sizes.
- Gauge size selection will depend on the type of medication to be infused and the flow rate required
- 19 gauge is the largest size PORT-A-CATH® needle offered. This gauge would be used for viscous solutions such as packed red blood cells or TPN.

Extension tubing is lipid-compatible.

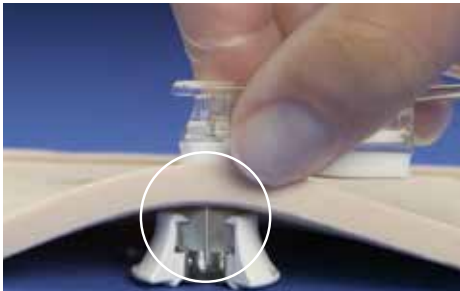


TOO SHORT

The needle is too short if it does not reach the reservoir floor of the portal. This may result in restriction of flow, obstruction of flow, or leakage of the fluid into the tissue surrounding the port.



TOO LONG

A GRIPPER PLUS® safety needle is too long for the portal if the safety arm does not rest flat on the base.

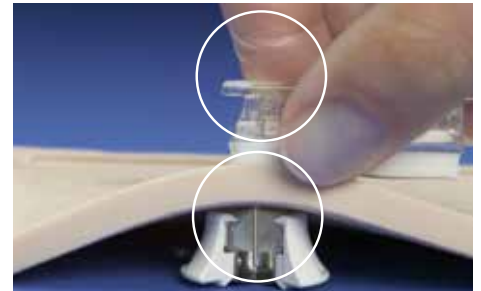

CORRECT LENGTH

The needle length is correct if the safety arm rests flat on the base and the needle touches the reservoir floor.

ORDERING INFORMATION

Needleless Y-site Luer-activated needleless injection site

	Length	Reorder No.	UPC No.	Unit
22 gauge	5/8 inch	21-2960-24	(01) 3 0610586 02768 4	12/box
	3/4 inch	21-2961-24	(01) 3 0610586 02769 1	12/box
	1 inch	21-2962-24	(01) 3 0610586 02770 7	12/box
	1 ¼ inch	21-2963-24	(01) 3 0610586 02771 4	12/box
20 gauge	5/8 inch	21-2964-24	(01) 3 0610586 02772 1	12/box
	3/4 inch	21-2965-24	(01) 3 0610586 02773 8	12/box
	1 inch	21-2966-24	(01) 3 0610586 02774 5	12/box
	1 ¼ inch	21-2967-24	(01) 3 0610586 02775 2	12/box
19 gauge	3/4 inch	21-2968-24	(01) 3 0610586 02776 9	12/box
	1 inch	21-2969-24	(01) 3 0610586 02777 6	12/box
	1 ¼ inch	21-2970-24	(01) 3 0610586 02778 3	12/box

Y-site Removable split septum injection cap on Luer-lock Y-site

	Length	Reorder No.	UPC No.	Unit
22 gauge	5/8 inch	21-2860-24	(01) 3 0610586 02564 2	12/box
	3/4 inch	21-2861-24	(01) 3 0610586 02565 9	12/box
	1 inch	21-2862-24	(01) 3 0610586 02566 6	12/box
	1 ¼ inch	21-2863-24	(01) 3 0610586 02567 3	12/box
20 gauge	5/8 inch	21-2864-24	(01) 3 0610586 02568 0	12/box
	3/4 inch	21-2865-24	(01) 3 0610586 02569 7	12/box
	1 inch	21-2866-24	(01) 3 0610586 02570 3	12/box
	1 ¼ inch	21-2867-24	(01) 3 0610586 02571 0	12/box
19 gauge	3/4 inch	21-2868-24	(01) 3 0610586 02572 7	12/box
	1 inch	21-2869-24	(01) 3 0610586 02573 4	12/box
	1 ¼ inch	21-2870-24	(01) 3 0610586 02574 1	12/box
Without Y-site	Length	Reorder No.	UPC No.	Unit
22 gauge	5/8 inch	21-2770-24	(01) 3 0610586 02563 5	12/box
	3/4 inch	21-2766-24	(01) 3 0610586 02559 8	12/box
	1 inch	21-2760-24	(01) 3 0610586 02553 6	12/box
	1 ¼ inch	21-2761-24	(01) 3 0610586 02554 3	12/box
20 gauge	5/8 inch	21-2769-24	(01) 3 0610586 02562 8	12/box
	3/4 inch	21-2767-24	(01) 3 0610586 02560 4	12/box
	1 inch	21-2762-24	(01) 3 0610586 02555 0	12/box
	1 ¼ inch	21-2763-24	(01) 3 0610586 02556 7	12/box
19 gauge	3/4 inch	21-2768-24	(01) 3 0610586 02561 1	12/box
	1 inch	21-2764-24	(01) 3 0610586 02557 4	12/box
	1 ¼ inch	21-2765-24	(01) 3 0610586 02558 1	12/box

PRODUCT(S) DESCRIBED MAY NOT BE LICENSED OR AVAILABLE FOR SALE IN CANADA AND OTHER COUNTRIES

Smiths Medical International
1500 Eureka Park, Lower Pemberton
Ashford Kent, TN25 4BF
Tel: +44 (0)845 850 0445

Smiths Medical ASD, Inc.
6000 Nathan Lane North
Minneapolis, MN 55442, USA
Tel: 1-614-210-7300
Toll-Free USA: 1-800-258-5361
www.smiths-medical.com

MPAUC-1410

Rx
ONLY
CE
2797

smiths medical

Find your local contact information at: www.smiths-medical.com/customer-support

Smiths Medical is part of the global technology business Smiths Group plc. Please see the Instructions for Use/Operator's Manual for a complete listing of the indications, contraindications, warnings and precautions. GRIPPER PLUS, GRIPPER, PORT-A-CATH, P.A.S. PORT, ProPort, and the Smiths Medical and Deltac design marks are trademarks of Smiths Medical. The symbol ® indicates the trademark is registered in the U.S. Patent and Trademark Office and certain other countries. All other names and marks mentioned are the trademarks or service marks of their respective owners. Product referenced is CE-marked. ©2019 Smiths Medical. All rights reserved. VA19374GB-102019