



APPEARANCE


Front View

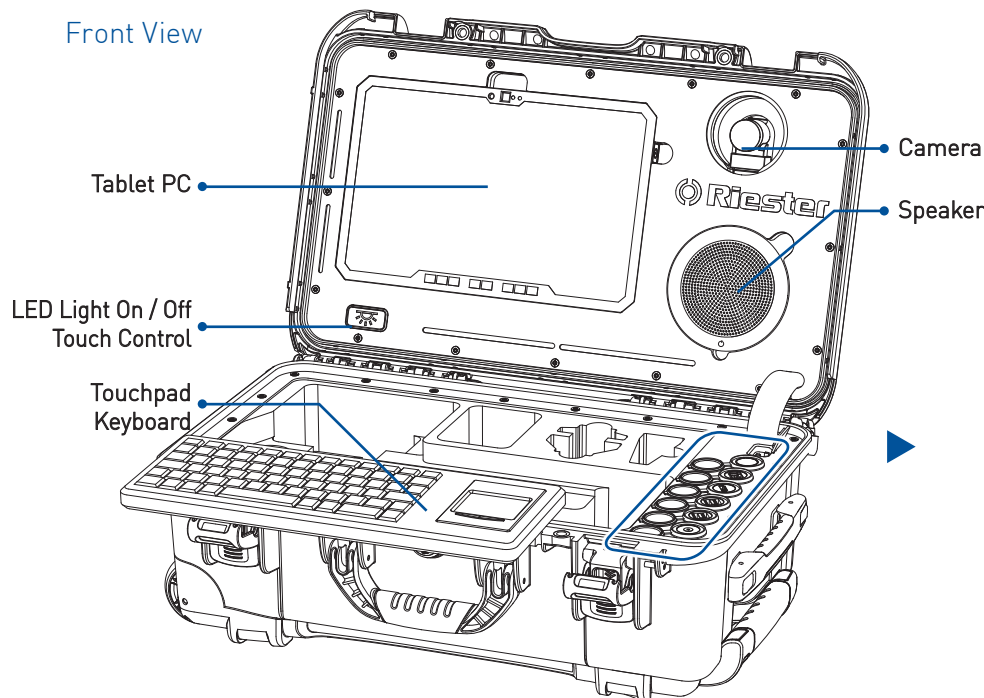

Detail View

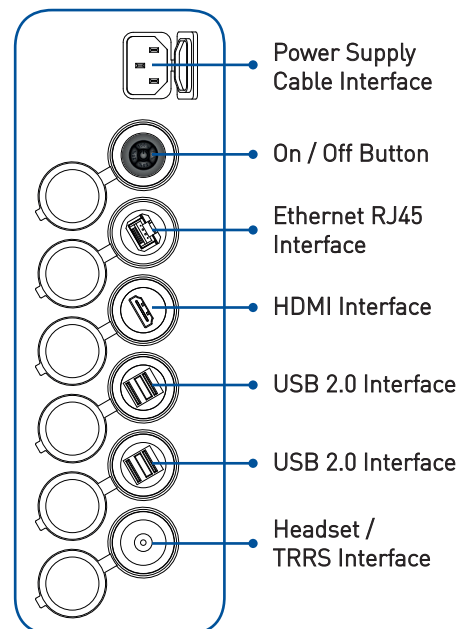

Medical Devices (integrated)

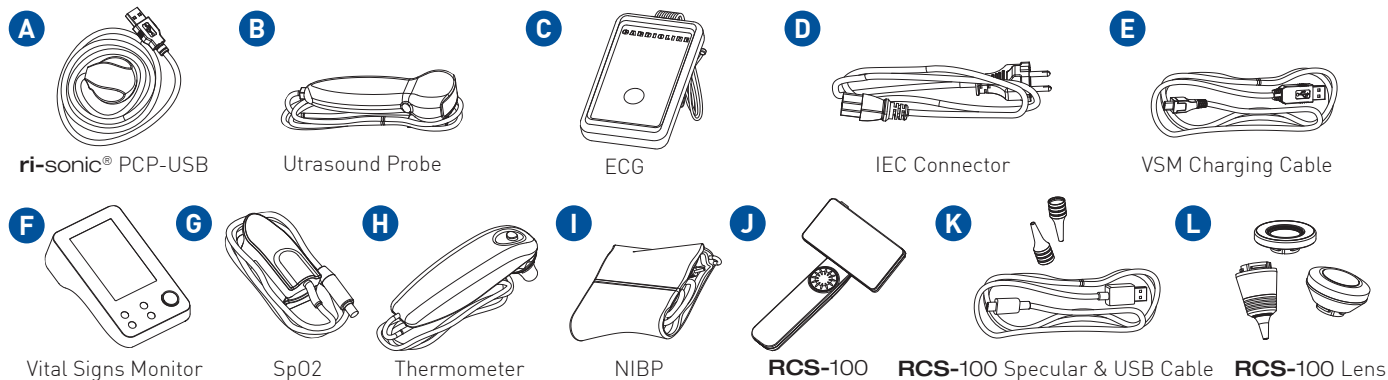

ORDERING INFORMATION


1920-EU	Telemedicine Case - EU	1920-UK	Telemedicine Case - UK
1920-AUS	Telemedicine Case - AUS	1920-US	Telemedicine Case - US
1920-HOSPITAL	Telemedicine Case - Hospital Grade US		
Medical Devices (Integrated)			
1840-BT	Non-contact thermometer ri-thermo® sensioPRO+ Bluetooth®	1966	PC-303 Portable Vital Signs Monitor
1926	Ultrasound Probe GP-C01	1970-HD60	RCS-100 System Medical Camera System
1927	touchECG System HD+12	4301	ri-sonic® PCP-USB Electronic Stethoscope
Medical Devices (Standalone)*			
13378	Spirometer Lung Monitor Bluetooth®	1923**	DEC200 Digital Ophthalmoscope Camera
1735-BT	Blood pressure device ri-champion® smartPRO+ Bluetooth®	1924**	LT-300 SD Digital Colposcope
1835-BT	Tympanic thermometer ri-thermo® tymPRO+ Bluetooth®	1928	Ultrasound Probe SP-L01

* In addition to the integrated devices, please note that you also have the option to order standalone devices. This means that besides the devices already included in the Telemedicine Case, you can place separate orders for additional devices based on your specific requirements.


**Riester's Telemedicine Case supports the add-on devices; however, it is important to note that these devices do not fit inside the case.

USING THE DEVICE

STEP 1. Open the Riester Telemedicine Case

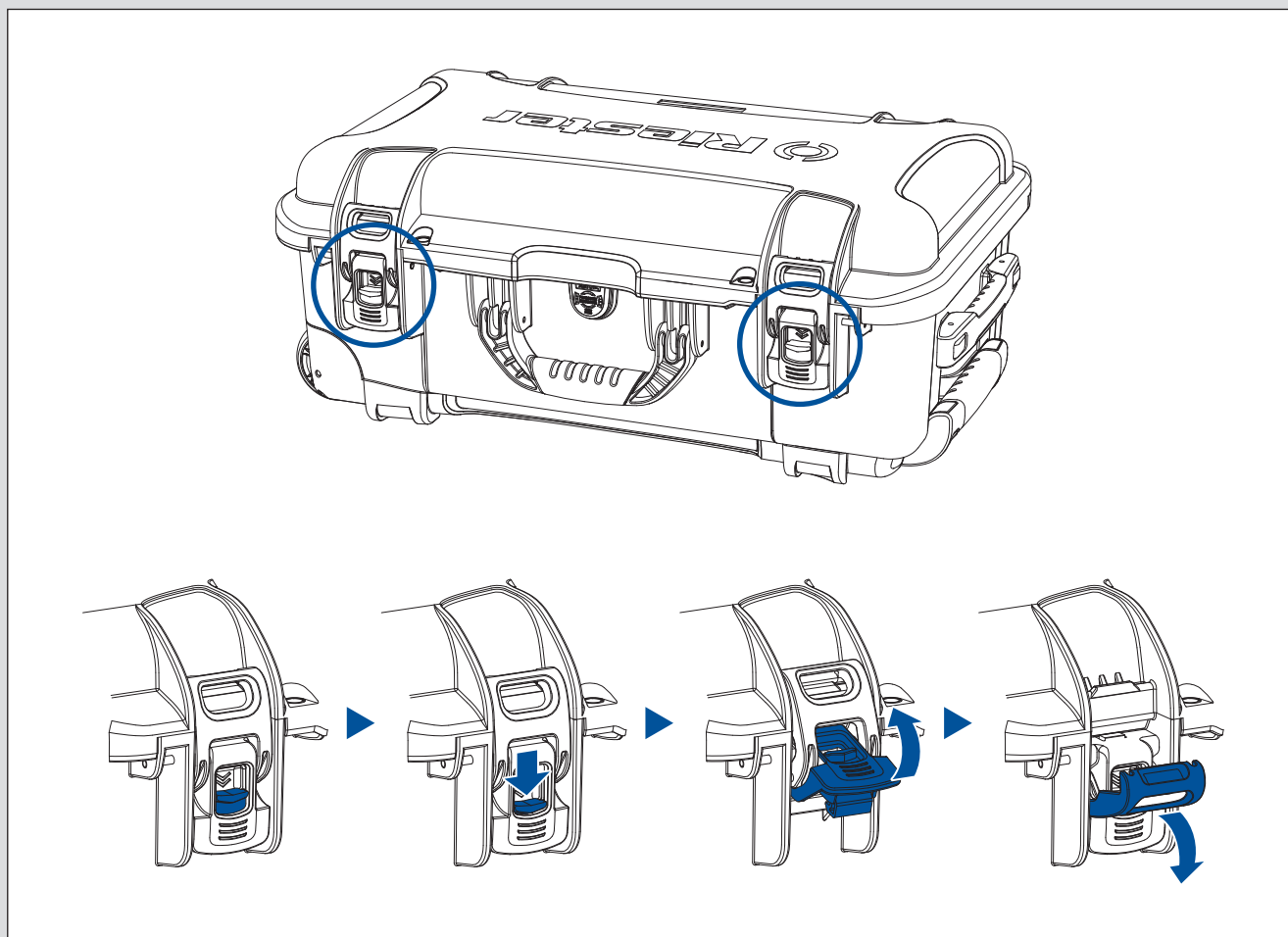


STEP 2. Telemedicine Case charging & Battery information

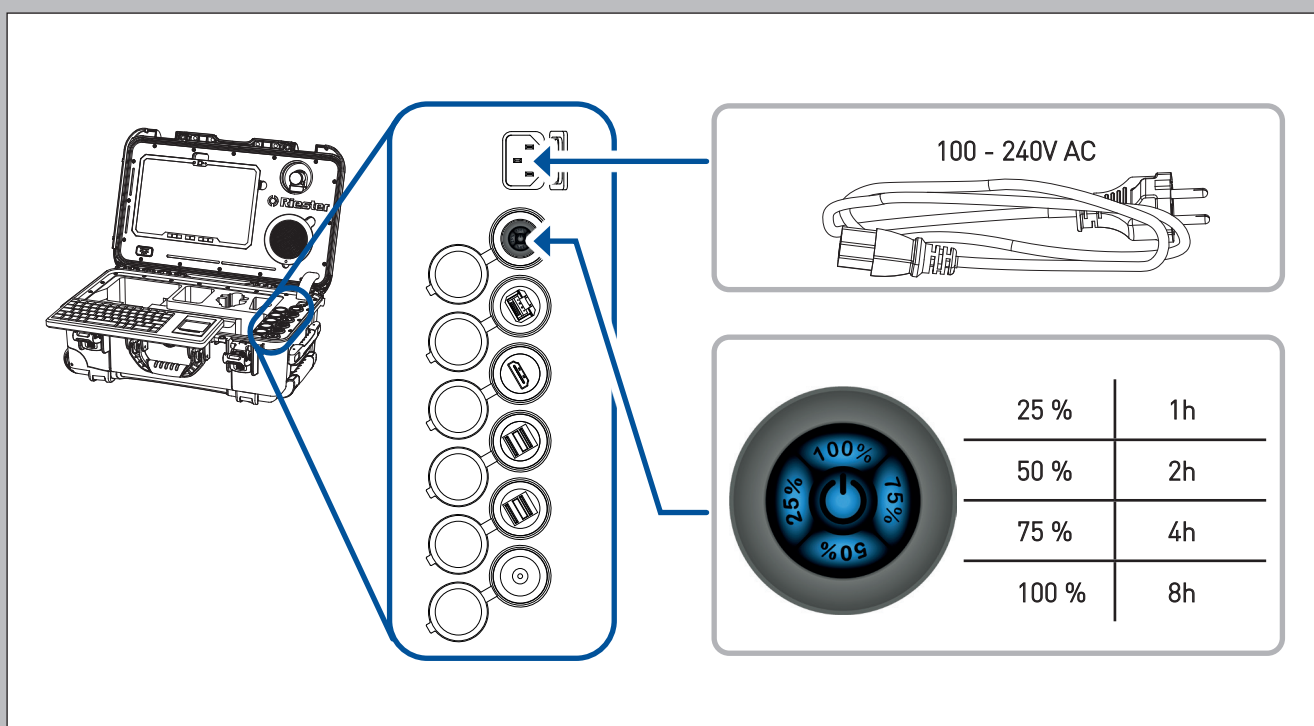


USING THE DEVICE

STEP 3. Physically place the relevant medical devices inside the Telemedicine Case

Bottom part of the Telemedicine Case

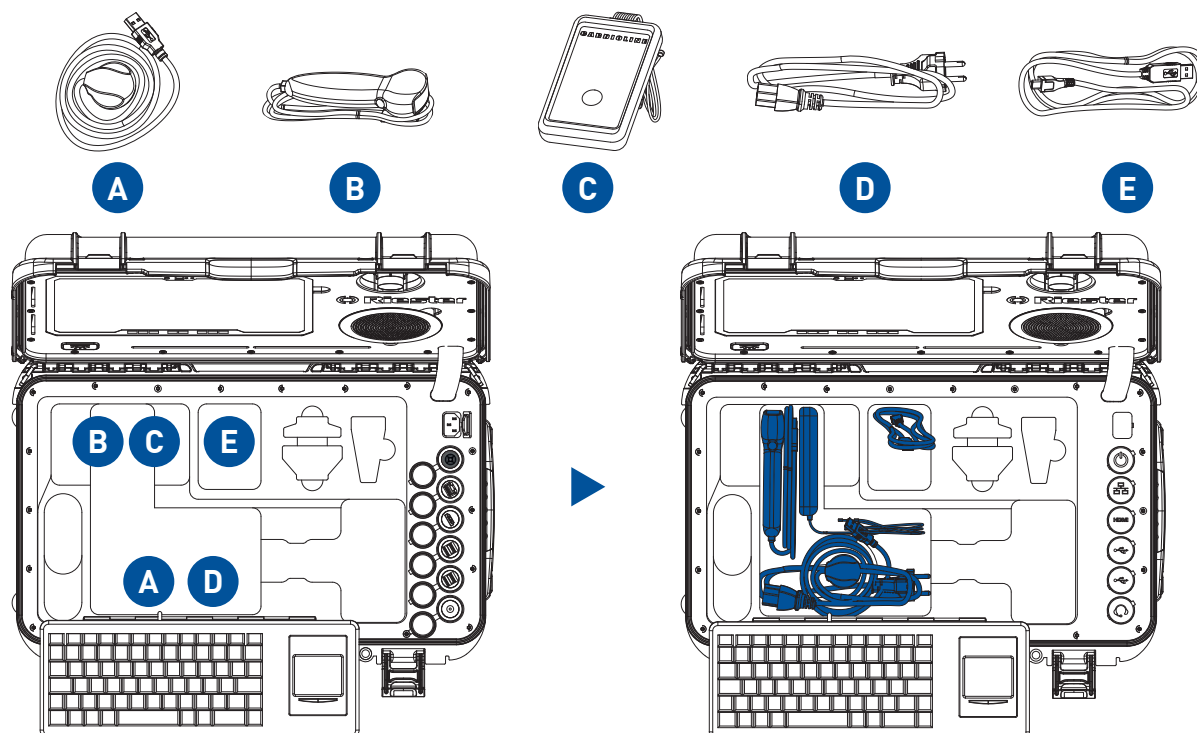


Top part of the Telemedicine Case

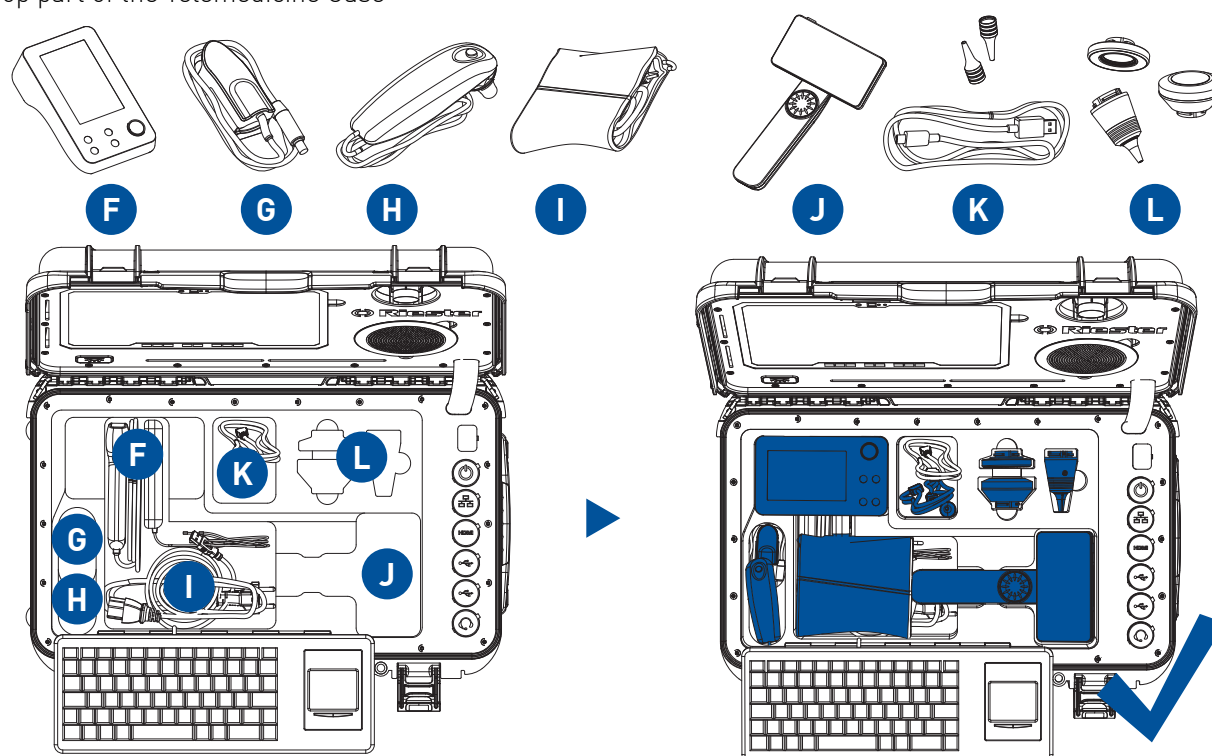


USING THE DEVICE

STEP 4. Turn on / off power

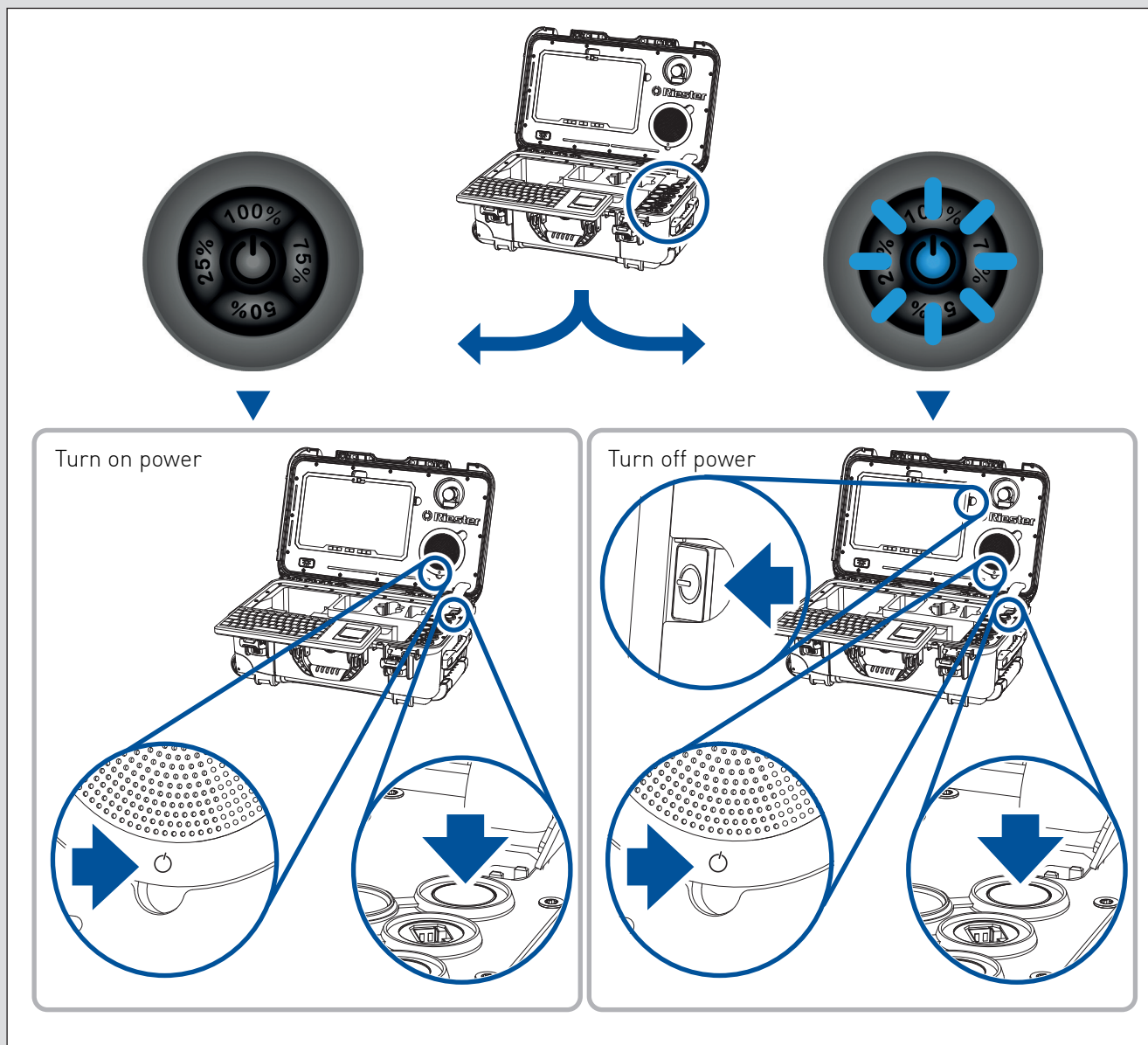


STEP 5. Check available Wi-Fi or connect ethernet

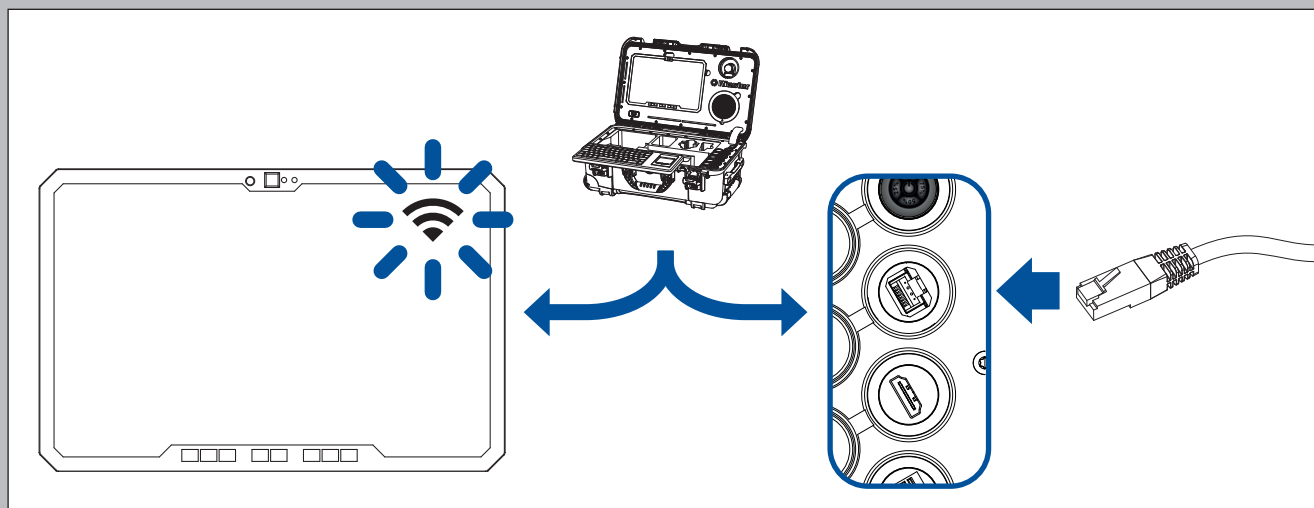


USING THE DEVICE

STEP 6. Connect headset / TRRS

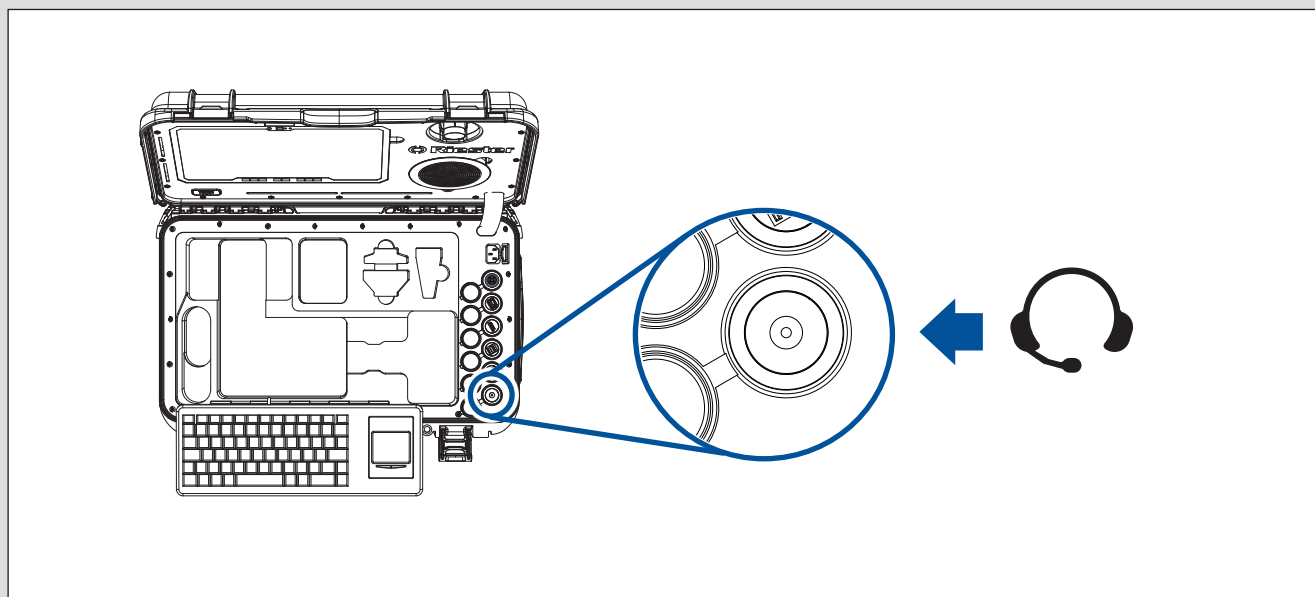


STEP 7. Connectivity medical devices

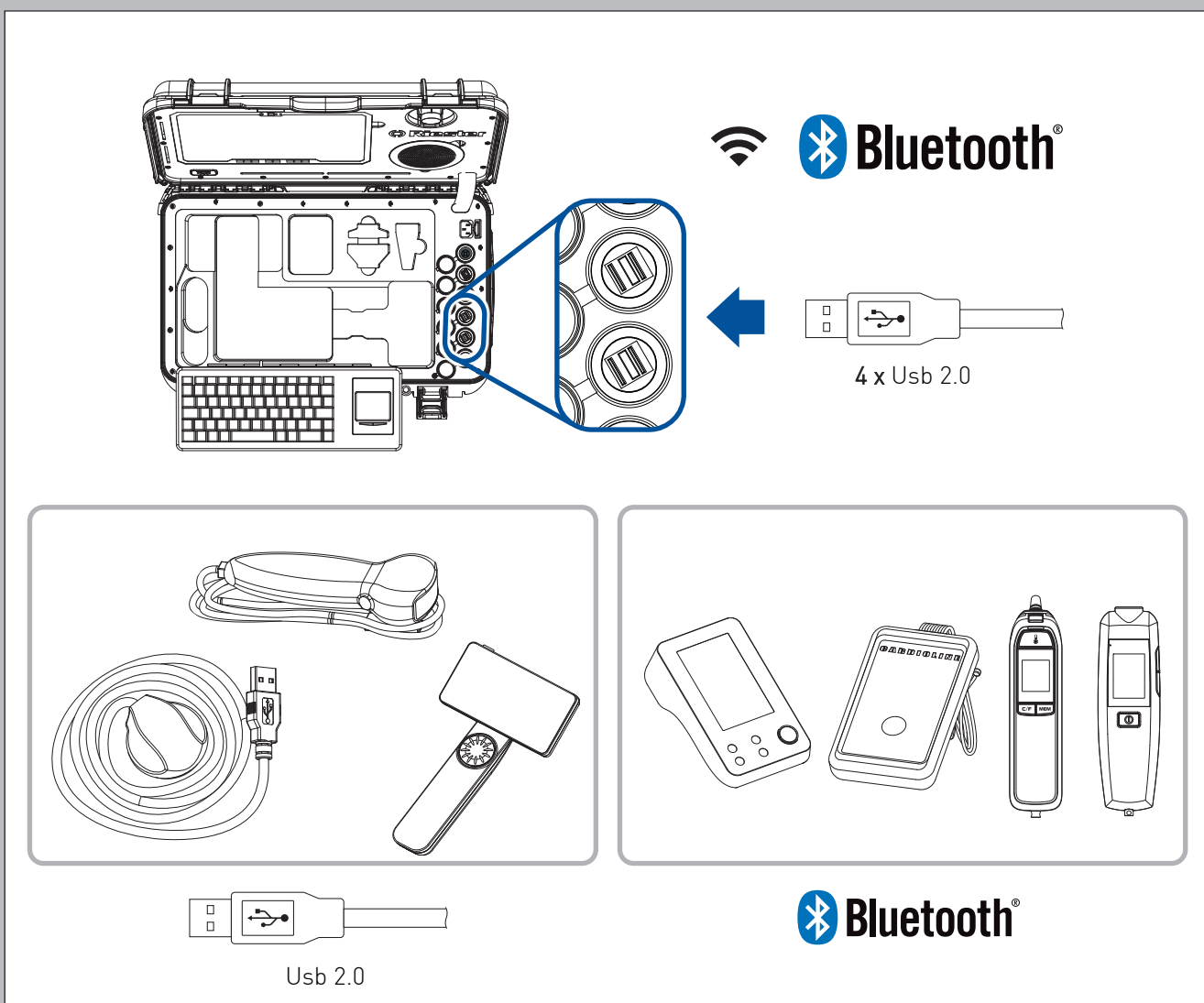


USING THE DEVICE

STEP 8. Connectivity Vital Signs Monitor – Connect the attached medical devices to the Vital Signs Monitor and connect the Vital Signs Monitor to the Telemedicine Case

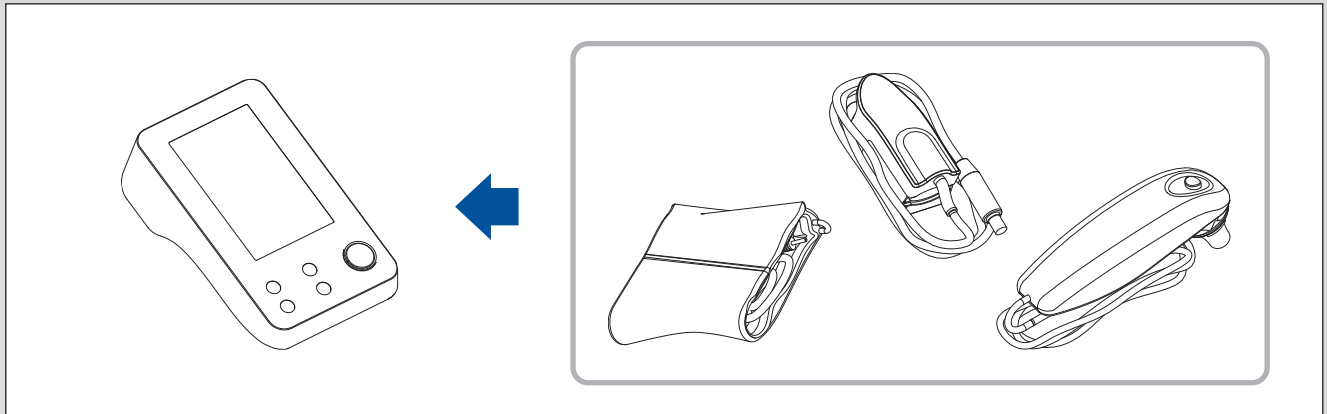

STEP 9. Log-in Teledoc

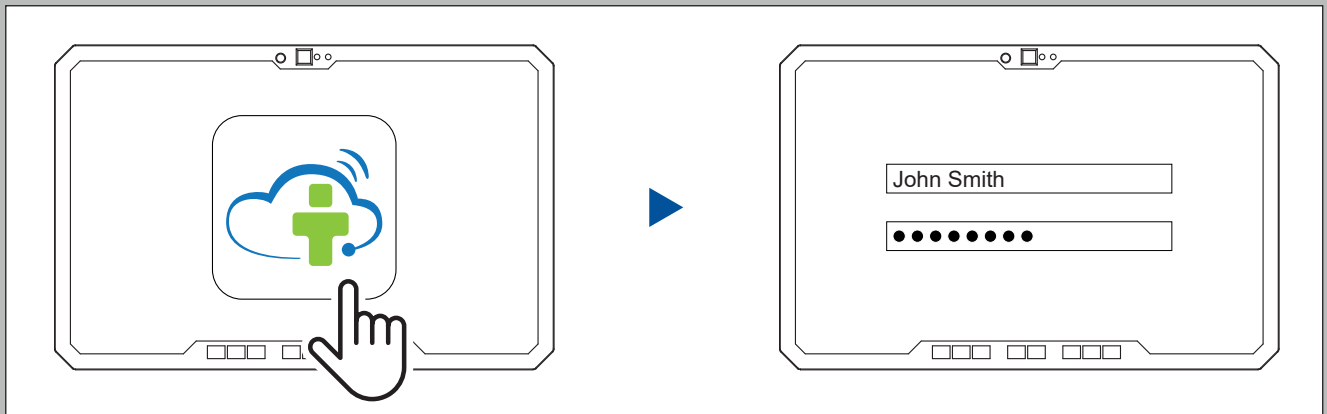

STEP 10. Further information



Please visit **riester.de** for more information



Refer to instructions for use for full functional usage



Rudolf Riester GmbH P.O. Box 35 | Bruckstraße 31 | 72417 Jungingen | Germany
Tel.: (+49) 7477-9270-0 | Fax.: (+49) 7477-9270-70 | info@riester.de | www.riester.de

Riester