



ASSUREWEAR™ CHEMO
TESTED VERSAGOWN®
WITH FLEXNECK™



TECHNICAL DATA SHEET



FEATURES & BENEFITS

- o VersaGown® with Flexneck™ technology significantly reduces cross-contamination versus standard gowns.
- o Designed using clinician feedback and independent clinical testing.
- o Provides nearly two-fold reduction in contamination of hand and wrist areas.
- o Helps ensure instant compliance—even without training.
- o Innovative design allows freedom of movement without compromising protection.
- o Flexneck™ creates the widest neck opening on the market and has ultrasonically sealed seams to prevent tearing during donning/doffing.
- o Winner of the Vizient Innovative Technology Award.

RECOMMENDED USE

- o Disposable Impervious Isolation Gown designed to protect health professionals and patients from cross infection.

CLAIM

- o Protective ultrasonically sealed gown that meets AAMI PB70 Level II barrier performance for protective appeal and tested for chemotherapy drug permeation as per ASTM D6978.

PRODUCT DETAILS

Product Name:	Ritmed® AssureWear™ Chemo Tested VersaGown® with Flexneck™ Technology
Product Type:	Single-use, non-sterile
Color:	See Appendix 1
Sizes:	See Appendix 1
Quantity per Bag:	10
Country of Origin:	China

DONNING INSTRUCTIONS

- o Perform hand hygiene before removing item from package.
- o Select appropriate garment for intended use (length, coverage, fluid resistance, etc.)
- o Fasten correctly: tied at neck and waist to ensure comfortable fit, adequate coverage and freedom of movement.
- o Protective apparel should not be worn in restrooms, personal offices or break rooms.

DOFFING INSTRUCTIONS

- o Remove apparel immediately after use. Do not re-use.
- o Perform proper hand hygiene after removing apparel to prevent cross-contamination.
- o Turn garment inside out when removing to keep infectious material inside.
- o Designate a safe receptacle in which to dispose of protective apparel after use.

CERTIFICATION & STANDARDS


ISO 13485 certified manufacturer.
Class 1 Medical Device.
Health Canada registered medical device manufacturer.

STORAGE CONDITIONS

These gowns should be disposed of after each use.
Store under ambient temperature and dry condition.



GOWN SPECIFICATIONS					
CHARACTERISTIC		SPECIFICATION			
Material Composition		Polypropylene (PP) & Polyethylene (PE). Latex-free			
Description		Gown made of non-woven material with full back closure. Features: inclined sleeves, flat belt, elastic cuffs (elastic does not go all around the cuff) with thumb loops available, Over-The-Head (OTH) flexible neck closure made of two elastic strings.			
Part of Gown		Dimension (cm)			
		Small	Medium	Large	X-Large
A	Gown Length	100	105	112	129
B	Underarm Width	63	67	71	79
	Open Gown Width	126	134	142	158
C	Underarm Height	63	68	76	90
D	Sleeve Length (Excludes Thumb Loop)	Up: 65 Down: 50	Up: 72 Down: 56	Up: 76 Down: 59	Up: 76 Down: 60
E	Shoulder Height	32	36	38	40
F	Shoulder Width	10	10	10	13
G1	Belt Length (Front)	87	87	87	87
G2	Belt Length (Back)	93	93	93	93
	Belt Width	5	5	5	5
	Belt Position	20 cm under the arm, 7 cm towards the front	20 cm under the arm, 7 cm towards the front	20 cm under the arm, 7 cm towards the front	20 cm under the arm, 7 cm towards the front
H	Cuff Width (Stretched)	14	14	14	14
I-1	Trim Length	54	54	54	54
I-2	Elastic Strings	9	9	9	9
I-3	Collar Height	10	10	10	10
I-4	Collar Width	19	19	22	22
	Tag Position	3,0 cm diameter, 2 cm from cuff edge	3,0 cm diameter, 2 cm from cuff edge	3,0 cm diameter, 2 cm from cuff edge	3,0 cm diameter, 2 cm from cuff edge
	Thumbloop Hole	23,5 (Left side of back only)	23,5 (Left side of back only)	23,5 (Left side of back only)	23,5 (Left side of back only)

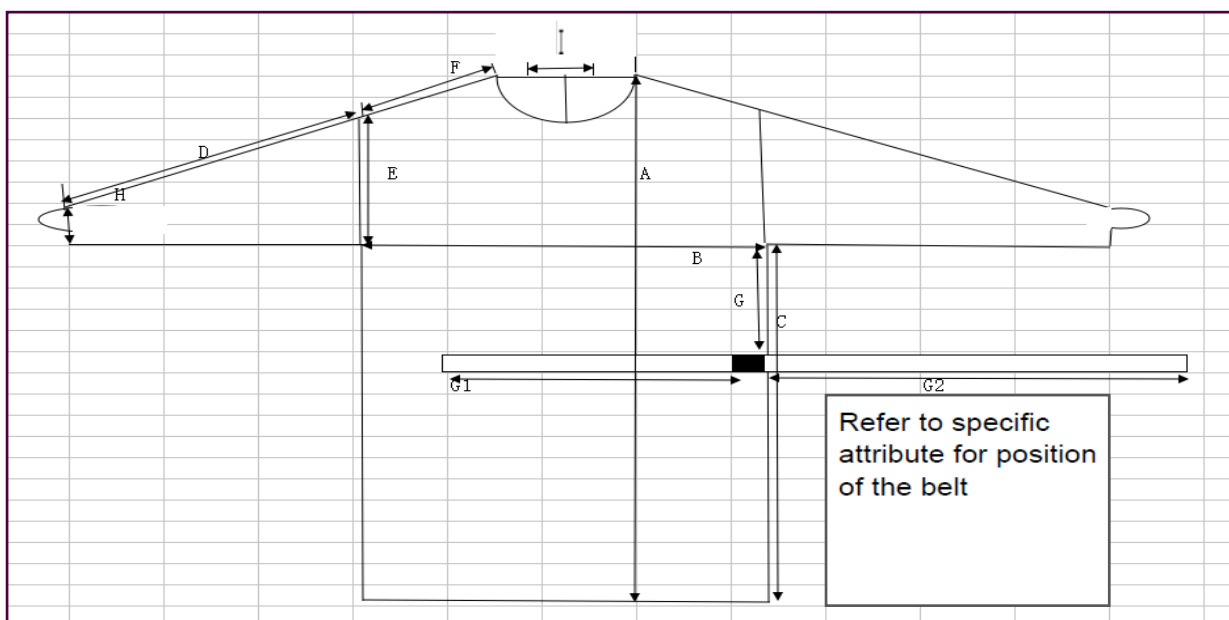

Appendix 1 (List of items)

Item Code	Description	Packaging Information	Size	Color	Specific Attribute	Treatment	Marketed in USA as	Marketed in Canada as
A69949	Ritmed AssureWear Chemo Tested VersaGown with Flexneck Technology	10 Gowns/Bag, 10 Bags/Case, 100 Gowns/Case	Small	Blue	Belt on the left side	Non-Sterile	Medical Device Class 1	Medical Device Class 1
A69950	Ritmed AssureWear Chemo Tested VersaGown with Flexneck Technology	10 Gowns/Bag, 10 Bags/Case, 100 Gowns/Case	Medium	Blue	Belt on the left side	Non-Sterile	Medical Device Class 1	Medical Device Class 1
A69951	Ritmed AssureWear Chemo Tested VersaGown with Flexneck Technology	10 Gowns/Bag, 10 Bags/Case, 100 Gowns/Case	Large	Yellow	Belt on the left side	Non-Sterile	Medical Device Class 1	Medical Device Class 1
A69954	Ritmed AssureWear Chemo Tested VersaGown with Flexneck Technology	10 Gowns/Bag, 10 Bags/Case, 100 Gowns/Case	Large	Blue	Belt on the left side	Non-Sterile	Medical Device Class 1	Medical Device Class 1
A69955	Ritmed AssureWear Chemo Tested VersaGown with Flexneck Technology	10 Gowns/Bag, 10 Bags/Case, 100 Gowns/Case	X-Large	Blue	Belt on the left side	Non-Sterile	Medical Device Class 1	Medical Device Class 1
A69954-C	Ritmed AssureWear Chemo Tested VersaGown with Flexneck Technology	10 Gowns/Bag, 10 Bags/Case, 100 Gowns/Case	Large	Blue	Belt on the left side	Non-Sterile	Medical Device Class 1	Medical Device Class 1

Appendix 1 (List of items)

Item Code	Description	Packaging Information	Size	Color	Specific Attribute	Treatment	Marketed in USA as	Marketed in Canada as
A69955-C	Ritmed AssureWear Chemo Tested VersaGown with Flexneck Technology	10 Gowns/Bag, 10 Bags/Case, 100 Gowns/Case	X-Large	Blue	Belt on the left side	Non-Sterile	Medical Device Class 1	Medical Device Class 1

AMD Medicom Inc. | 2555 Chemin de l'Aviation, Montreal, Quebec, Canada H9P 2Z2